

REPORT DOCUMENTATION PAGE

AFRL-SR-AR-TR-03-

Public reporting burden for this collection of information is estimated to average 1 hour per response, including the time for reviewing the data needed, and completing and reviewing this collection of information. Send comments regarding this burden estimate or at reducing this burden to Department of Defense, Washington Headquarters Services, Directorate for Information Operations and Research, 22202-4302. Respondents should be aware that notwithstanding any other provision of law, no person shall be subject to any penalty for not providing information if it does not display a currently valid OMB control number. PLEASE DO NOT RETURN YOUR FORM TO THE ABOVE ADDRESS.

1. REPORT DATE (DD-MM-YYYY) 2/26/03		2. REPORT TYPE Final Performance Report		3. DATES COVERED (From - To) 6/01/01 - 11/30/02	
4. TITLE AND SUBTITLE Nanotechnology Instrumentation				5a. CONTRACT NUMBER	
				5b. GRANT NUMBER F496200110419	
				5c. PROGRAM ELEMENT NUMBER	
				5d. PROJECT NUMBER	
6. AUTHOR(S) H. M. Gibbs, Ph.D and G. Khitrova, Ph.D.				5e. TASK NUMBER	
				5f. WORK UNIT NUMBER	
				8. PERFORMING ORGANIZATION REPORT NUMBER	
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) Optical Sciences Center The University of Arizona 1630 E. University Blvd. Tucson, Arizona 85721				10. SPONSOR/MONITOR'S ACRONYM(S)	
				11. SPONSOR/MONITOR'S REPORT NUMBER(S)	
9. SPONSORING / MONITORING AGENCY NAME(S) AND ADDRESS(ES) AFOSR/NE Dr. Gernot Pomrenke 4015 Wilson Boulevard, Room 713 Arlington, VA 22203-1954					
12. DISTRIBUTION / AVAILABILITY STATEMENT Approved for Public release - Distribution limited					
13. SUPPLEMENTARY NOTES N/A					
14. ABSTRACT Various instruments have been purchased and tested that extend existing nanotechnology capabilities and complement existing equipment. The fact that these laboratories have worked for several years on 3D microcavities and micron-size spectroscopic measurements clarified the need for several complementary instruments. The major new thrust made possible by this grant is spectroscopic capability at 1300 nm. A CCD camera system provides very sensitive detection from 800 to 1600 nm. The existing fs Ti:Sa laser was modified to optimize it for pumping a new optical parametric oscillator, providing short pulses from 1100 to 1600 nm, and the new Millennium X solid state pump provides improved beam stability and efficiency. These instruments have enabled the study of quantum dots and 3D nanocavities, both photonic-crystals and microdisks, in the 1000 - 1300 nm range. A cryostat with nanopositioners within the vacuum has greatly facilitated these measurements, requiring micron stability for minutes. Its temperature control enables scanning of the dot-nanocavity detuning and temperature-dependent measurements.					
15. SUBJECT TERMS Quantum dots, nanocavities emitting at about 1300 nm at room temperature.					
16. SECURITY CLASSIFICATION OF:			17. LIMITATION OF ABSTRACT	18. NUMBER OF PAGES 9	19a. NAME OF RESPONSIBLE PERSON Prof. Hyatt Gibbs
a. REPORT	b. ABSTRACT	c. THIS PAGE			19b. TELEPHONE NUMBER (include area code) (520) 621-2941

20030515 113

Nanotechnology Instrumentation

Final Performance Report February 2003

Hyatt M. Gibbs, Ph.D,
Optical Sciences Center
University of Arizona
Tucson, Arizona 85721
gibbs@optics.arizona.edu

Galina Khitrova, Ph.D
Optical Sciences Center
University of Arizona
Tucson, Arizona 85721
galina@optics.arizona.edu

**US Air Force Office of Scientific Research
Grant No. F49620-01-1-0419**

**APPROVED FOR PUBLIC RELEASE
DISTRIBUTION LIMITED**

THE VIEWS, OPINIONS, AND/OR FINDINGS CONTAINED IN THIS REPORT ARE THOSE OF THE AUTHOR (S)
AND SHOULD NOT BE CONSTRUED AS AN OFFICIAL POSITION, POLICY
OR DECISION, UNLESS SO DESIGNATED BY OTHER DOCUMENTATION,

4. Accomplishments/New Findings

The actual purchases made are the following. The major items are as proposed. But some savings and slight modifications permitted additional cryostat purchases.

Vendor	Description	PO#	Amount
Roper Scientific	Princeton Instruments Digital CCD Camera System	P655182	\$ 30,690.45
Spectra Physics	High Power Diode Pumped Solid State Laser, etc	P610361	\$166,000.00
Janis Research	³ He Cryostat	P618898	\$ 11,600.00
Varian, Inc	Triscoll 300 1 Phase Current Pump	P615539	\$ 4,701.22
Hamamatsu	Assembly GaAsP Photon Counting Head with Cooler	P622643	\$ 983.00
CryoVac	Konti-Cryostat for Microscopic Measurements	P615431	\$ 26,212.80
TOTAL			\$240,187.47

The primary objective was to acquire new instruments enabling us to study single quantum dots and 3D nanocavities in the 1000 -1300 nm wavelength range (Fig. 1).

Objective	Study single quantum dots and 3D nanocavities. Couple quantum dots to a confined light field, demonstrate true strong coupling and quantum entanglement
Approach	Control and interrogate quantum dots and quantum-dot nanocavities as e.g. a strongly confined InAs quantum dot in a photonic-crystal point-defect nanocavity.
Relevance	The reversible regime of quantum mechanics can be useful for quantum information exchange and can lead to novel high-speed optoelectronic devices.

Fig. 1: Objective, approach and relevance.

In order of descending expenses, the instruments acquired were: Millennium X pump laser for the fs Ti:Sa laser, OPAL optical parametric oscillator, and corresponding modification of the fs Ti:Sa laser; 800 – 1600 nm InGaAs linear detector array; cryostat with internal nanopositioners, 16% purchase of ³He 0.4K cryostat; vacuum pump for pumping cryostats; and single-photon-counting photomultiplier. See Fig 2.

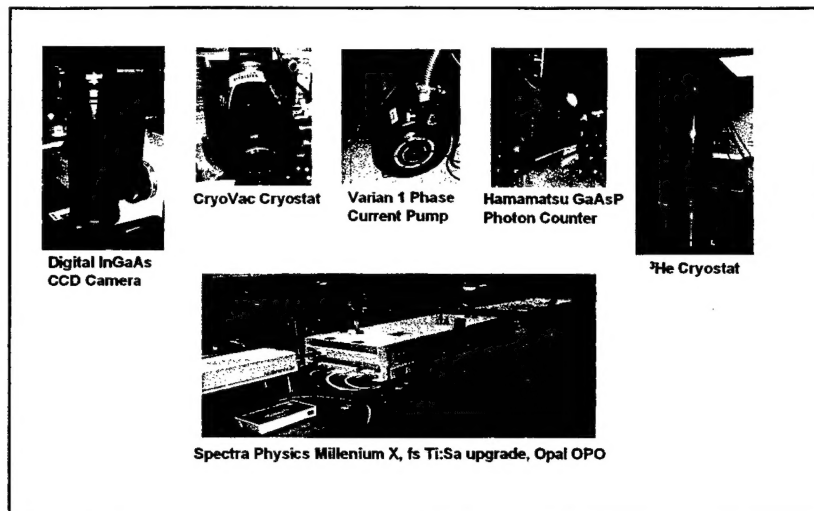


Fig. 2: Photographs of the purchased equipment.

The InGaAs detector array has been used to detect photoluminescence (PL) from individual quantum dots and photonic crystal nanocavities (Fig. 3) as well as 300K and 8K lasing of a quantum-dot microdisk (Fig. 4). The CryoVac cryostat enables its case and helium transfer lines to be secured to the table while the sample can be translated in two dimensions by internal nanopositioners (Fig. 5). It has been used for the PL (Fig. 3) and lasing (Fig. 4) experiments, for surveying large areas of samples in search of photonic-crystal cavities (Fig. 6, left), and for measuring the increase in quantum-well linewidth with temperature (Fig. 6, right).

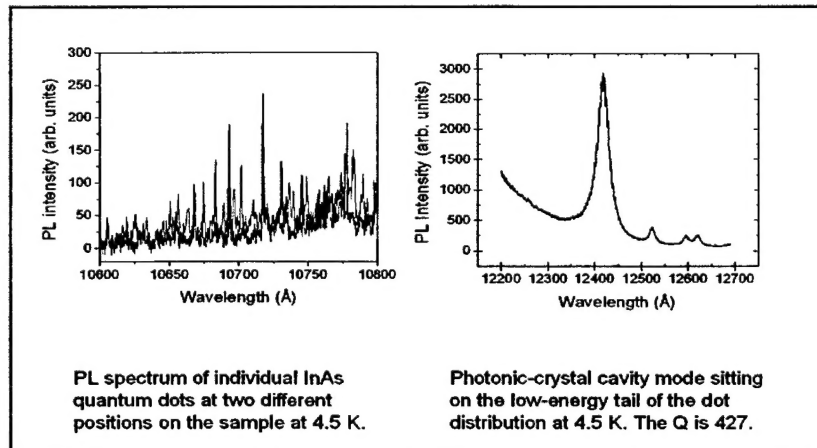


Fig. 3: PL from individual quantum dots and photonic crystal nanocavities.

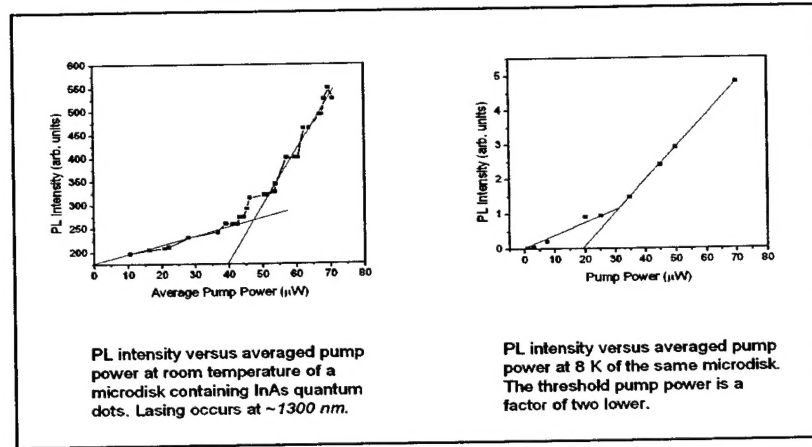


Fig. 4: Room temperature and 8K lasing of a quantum-dot microdisk.

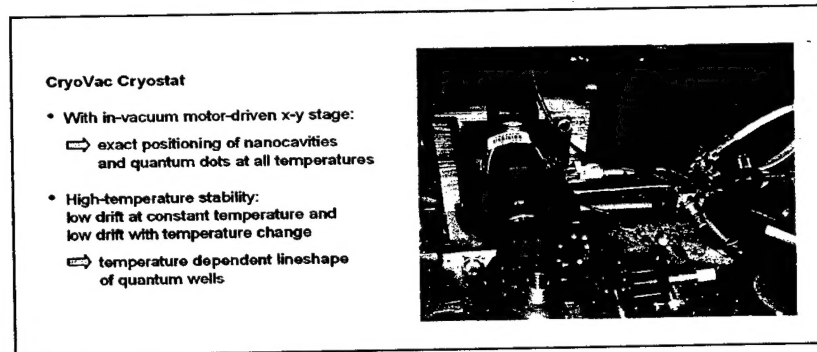


Fig. 5: CryoVac cryostat with internal x-y nanopositioners.


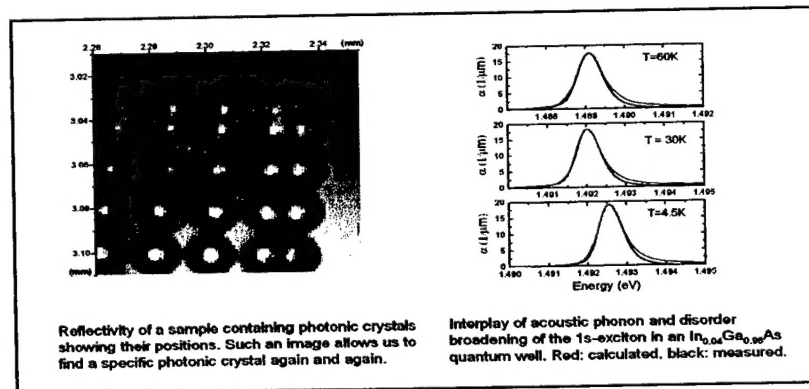


Fig. 6: 2D image of photonic crystals (left), and temperature-dependent change of the excitonic lineshape of an InGaAs quantum well (right).

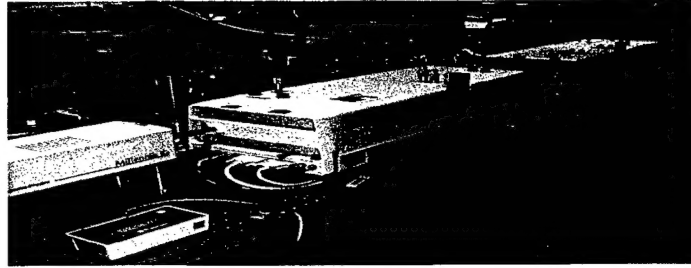

Spectra Physics Millennium X, fs Ti:Sa upgrade, Opal OPO

Fig. 7: Femtosecond laser system: solid-state pump, 100fs Ti:Sa oscillator and opto-parametric oscillator (OPO).

The Millennium X is an all-solid-state laser providing up to 10 W at 530 nm to pump the fs Ti:Sa oscillator (Fig. 7). It has far superior pointing stability, much higher wall-plug efficiency, and lower maintenance cost than the ten-year-old argon laser it replaced. This pump-oscillator pair has been used for preliminary studies on PL from InGaAs quantum wells (Fig. 8). This study has delayed the use of the OPAL OPO to excite quantum dot sample for lifetime measurements by upconversion.

Objective	Determine what fraction of carriers form excitons before recombining in InGaAs QWs after nonresonant excitation.
Motivation	PL at 1s exciton resonance maybe mostly electron-hole plasma emission, rather than excitonic.
Approach	Ratio $PL(E)/absorption(E)$ allows to determine deviation from purely plasma PL.
Relevance	Possible exciton condensation, when excitonic population is large enough and lifetime is long enough.

Fig. 8: Studies on PL from InGaAs quantum wells.

The Hamamatsu photon-counting photomultiplier has been used to determine the lifetimes of type-II excitons in a GaAs/AlAs superlattice structure (Fig. 9). By obtaining manufacturer discounts and eliminating the OPAL doubler that we can build ourselves, we had \$11,600.00 left that we added to other funds to purchase a ^3He cryostat (Fig. 10). Before we could reach just below 2K by pumping on ^4He ; the ^3He cryostat can go below 400 mK. This new capability may be crucial if exciton condensation is to be achieved.

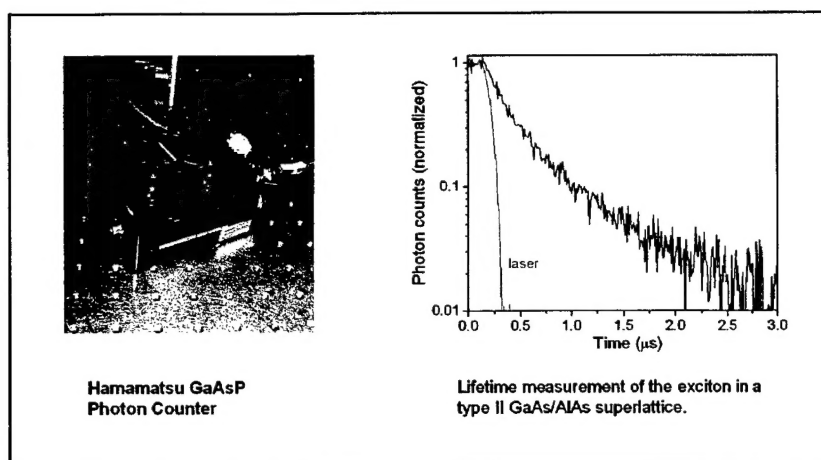


Fig. 9: Long-lifetime measurement of an exciton in a type II superlattice.

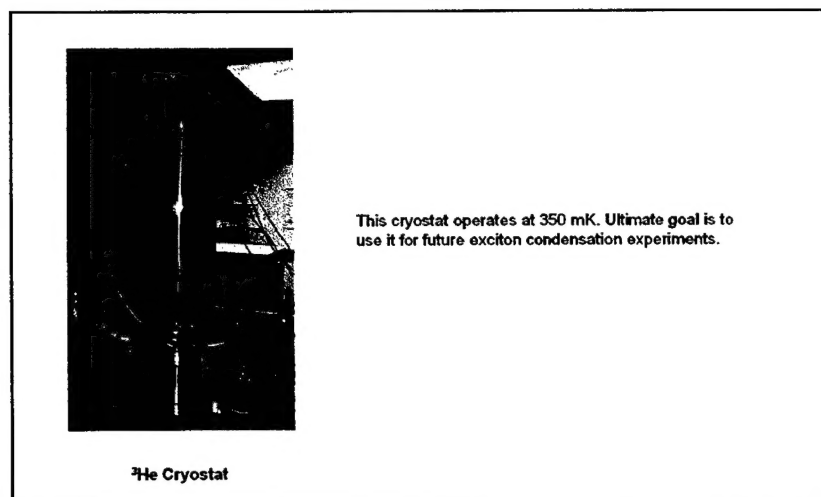


Fig. 10: ^3He cryostat for optical measurements at temperatures below 400 mK.

Experiments on single quantum dots and submicron-diameter photonic crystal nanocavities are tedious and time consuming even with these new state-of-the-art instruments. This is because of the difficulty of isolating the samples from environmental perturbations that can move them by more than a micron in one to five minutes, typical integration times for these weak signals. Nonetheless we are succeeding in seeing PL from single quantum dots and 3D nanocavities. All of our runs at present center on studying InAs quantum dot samples grown by Prof. Dennis Deppe, University Texas at Austin, and fabricated into phonic-crystal cavities by Prof. Axel Scherer and Tomo Yoshie, of Caltech, or microdisk cavities by Prof. John O'Brien, USC. The primary goal

is to see strong coupling between a single quantum dot and a single cavity mode, characterized by a double-peaked PL spectrum that should exhibit an anti-crossing as the quantum dot transition is temperature scanned through the cavity resonance. So far the cavity linewidth has been larger than the calculated splitting, or the dot density has been too high to isolate a single dot. A sample with a single layer of dots has been grown and will soon be processed and sent to us. The instruments purchased with this grant are essential to this research program.

5. Personnel supported

None supported directly by this equipment grant, but associated with the use of the instruments are:

Professors H. M. Gibbs and G. Khitrova
Research Professors C. Ell and J. Xu
Graduate Students G. Rupper, S. Chatterjee, and S. Mosor

6. Publications

A. Thränhardt, C. Ell, S. Mosor, G. Rupper, G. Khitrova, H.M. Gibbs, and S. W. Koch, "Interplay of phonon and disorder scattering in semiconductor quantum wells", submitted to Physical Review B, 2003.

7. Interaction/Transitions

W. Hoyer, M. Kira, S. W. Koch, P. Brick, S. Chatterjee, C. Ell, G. Khitrova, and H. M. Gibbs, "Nonequilibrium characteristics of excitonic luminescence", in *OSA Trends in Optics and Photonics (TOPS) Vol. 74, Quantum Electronics and Laser Science Conference (QELS 2002)*, Technical Digest, Postconference Edition (Optical Society of America, Washington DC, 2002, p. 106).

S. Chatterjee, C. Ell, G. Khitrova, H. M. Gibbs, W. Hoyer, M. Kira, and S. W. Koch, "Exciton formation in semiconductor quantum wells", Seminar talk at the University of Marburg, Germany, 2002.

S. Chatterjee, S. Mosor, C. Ell, G. Khitrova, H. M. Gibbs, W. Hoyer, M. Kira, and S. W. Koch, "Exciton formation in semiconductor quantum wells", Poster at the Photonics Initiative Workshop, Tucson, 2003.

8. New discoveries, inventions, or patent disclosures

None.

9. Honors/Awards:

H. M. Gibbs: Michelson Medal 1994.

H. M. Gibbs: Humboldt Research Award 1998.